

# Targeted Proteomics coming of age – SRM PRM and DIA performance evaluated from a core facility perspective

Tobias Kockmann<sup>1</sup>, Christian Trachsel<sup>2</sup>, Christian Panse<sup>1</sup>, Åsa Wähländer<sup>1</sup>, Nathalie Selevsek<sup>1</sup>, Jonas Grossmann<sup>1</sup>, Witold E. Wolski<sup>1</sup>, Claudia Fortes<sup>1</sup>, Paolo Nanni<sup>1</sup> and Ralph Schlapbach<sup>1</sup>

## 1 Introduction

Quantitative mass spectrometry is a rapidly evolving methodology applied in a large number of -omics type research projects. During the past years, new designs of mass spectrometers have been developed and launched as commercial systems while in parallel new data acquisition schemes and data analysis paradigms have been introduced. Core facilities provide access to such technologies, but also actively support the researchers in finding and applying the best-suited analytical approach. In order to implement a solid fundament for this decision making process, facilities need to constantly compare and benchmark the various approaches. In this work we compare the quantitative accuracy and precision of current label-free targeted proteomics approaches (SRM, PRM and DIA) across multiple liquid chromatography mass spectrometry (LC-MS) platforms, using a readily available commercial standard sample.

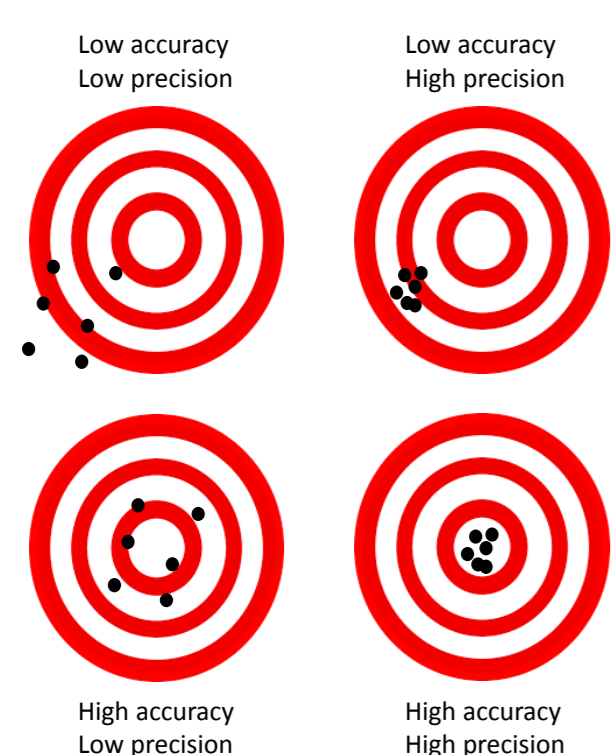


Figure 1: Accuracy and precision

## 2 Methods

A pre-digested protein sample (MSQC1; Sigma-Aldrich) with *a priori* known analyte quantities (mixture of 6 trypsin-digested human proteins and 14 corresponding heavy synthetic peptides) was targeted in a complex yeast matrix on all platforms. In a first part, the quantitative accuracy and precision (see Figure 1) was investigated at constant analyte quantities. In a second part, the impact of decreasing analyte quantities on analytical performance was studied over a concentration range of three orders of magnitude (dilution series data). Finally, the human influence on quantitative accuracy was investigated and contrasted with machine learning solutions (user study).

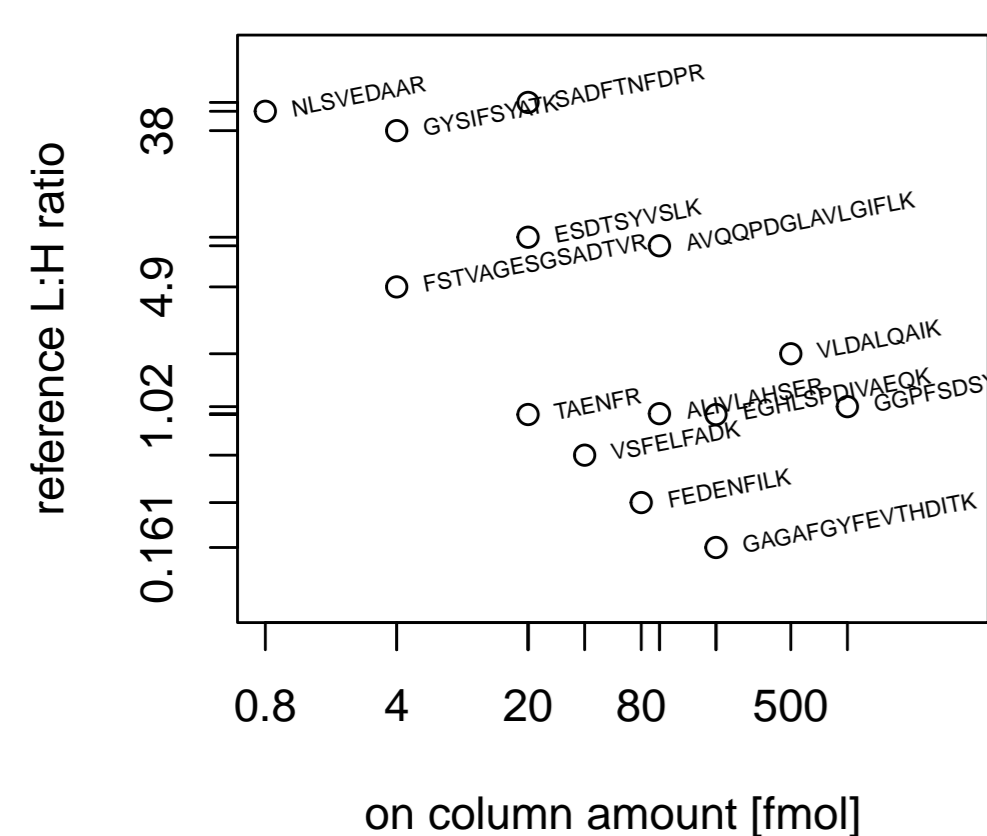


Figure 2: Target peptides – The scatter-plot displays the reference L:H ratio versus the on-column amount of heavy peptide. Note, x and y axis are drawn in log scale.

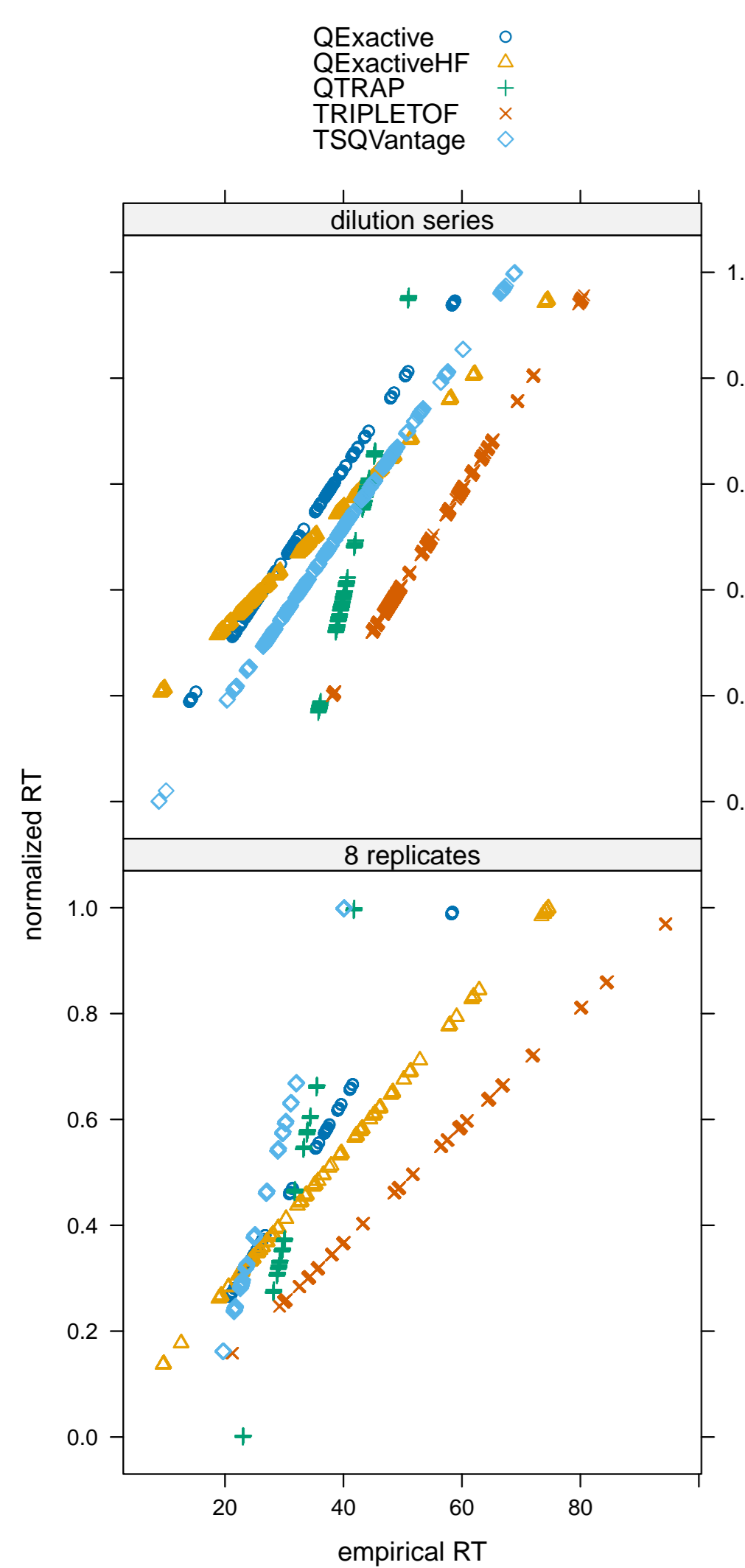


Figure 3: Chromatography – The graphs compare the LC gradient of each platform by plotting the normalized RT values versus the empirical RT values for the 8 replicate (bottom) and dilution series (top) data.

MS system	Analysis mode	PCS	IT	MS2 res.
TSQ Vantage	SRM	0.7 Da	20 ms	-
QTRAP 5500	SRM	0.7 Da	20 ms	-
Q EXACTIVE	PRM	2.0 Da	120 ms <sup>a</sup>	70'000 <sup>b</sup>
Q EXACTIVE HF	DIA	25.0 Da	45 ms <sup>c</sup>	30'000 <sup>d</sup>
TripleTOF 5600	SWATH	25.0 Da	100 ms	15'000

<sup>a</sup>or 1e5

<sup>b</sup>at 200 m/z

<sup>c</sup>or 3e6

<sup>d</sup>at 200 m/z

Table 1: **Measurement schemes** – Summarizes the measurement schemes compared in this study. PCS: Precursor selectivity, IT: Injection time, MS2 res.: resolution of MS2 scan. For Orbitrap type analyzers the max. injection time and the automated gain control (AGC) value is given. Filling of the C-trap will end as soon as one of the two parameters (IT or AGC) is fulfilled.

## 3 Results

– Good linearity for all targeted peptides across the tested concentration range was observed.

– Targeted analysis platforms delivered very reproducible (ie. intra-assay CVs less than 15%) and consistent quantitative data sets, even between platforms (see Figure 5).

– Expected fold changes were obtained for all dilution points, however at low peptide concentrations, the accuracy of the measurements were compromised to different extents depending on the MS platform (see Figure 6).

– Compared to the automated evaluation of peak-groups in Skyline a manual curation usually leads to a higher number of quantified peaks but only expert user manage to do this without compromising on the variance (see Figure 4). Each manual peakgroup validation will introduce a human bias which is not reproducible and therefore automated validation is favored (see Figure 4).

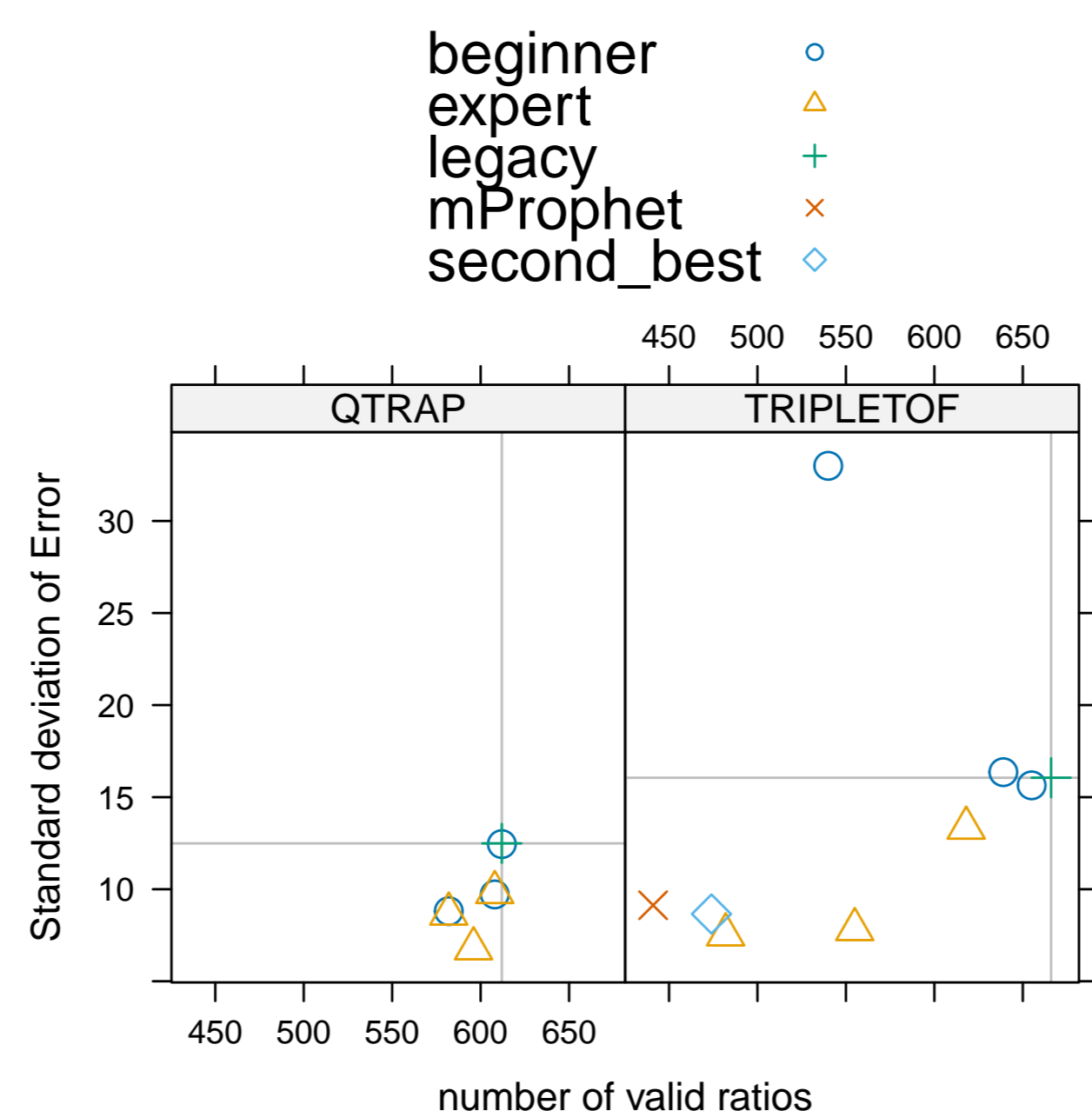


Figure 4: **Human impact on quantitative accuracy** – The scatterplots graph the standard deviation of the error between measured and reference  $\log_2$  L:H ratio against the number of valid ratios (L and H value is not NA) for each subject. The crossing grey lines indicate the algorithmic proposed start solution of the Skyline legacy peak picking. On the Q-trap no automated peak group selection was trained since no decoy transitions were measured. mProphet and second-best were trained on the 8 repl. data set and applied with a q-value cut-off of 0.001.


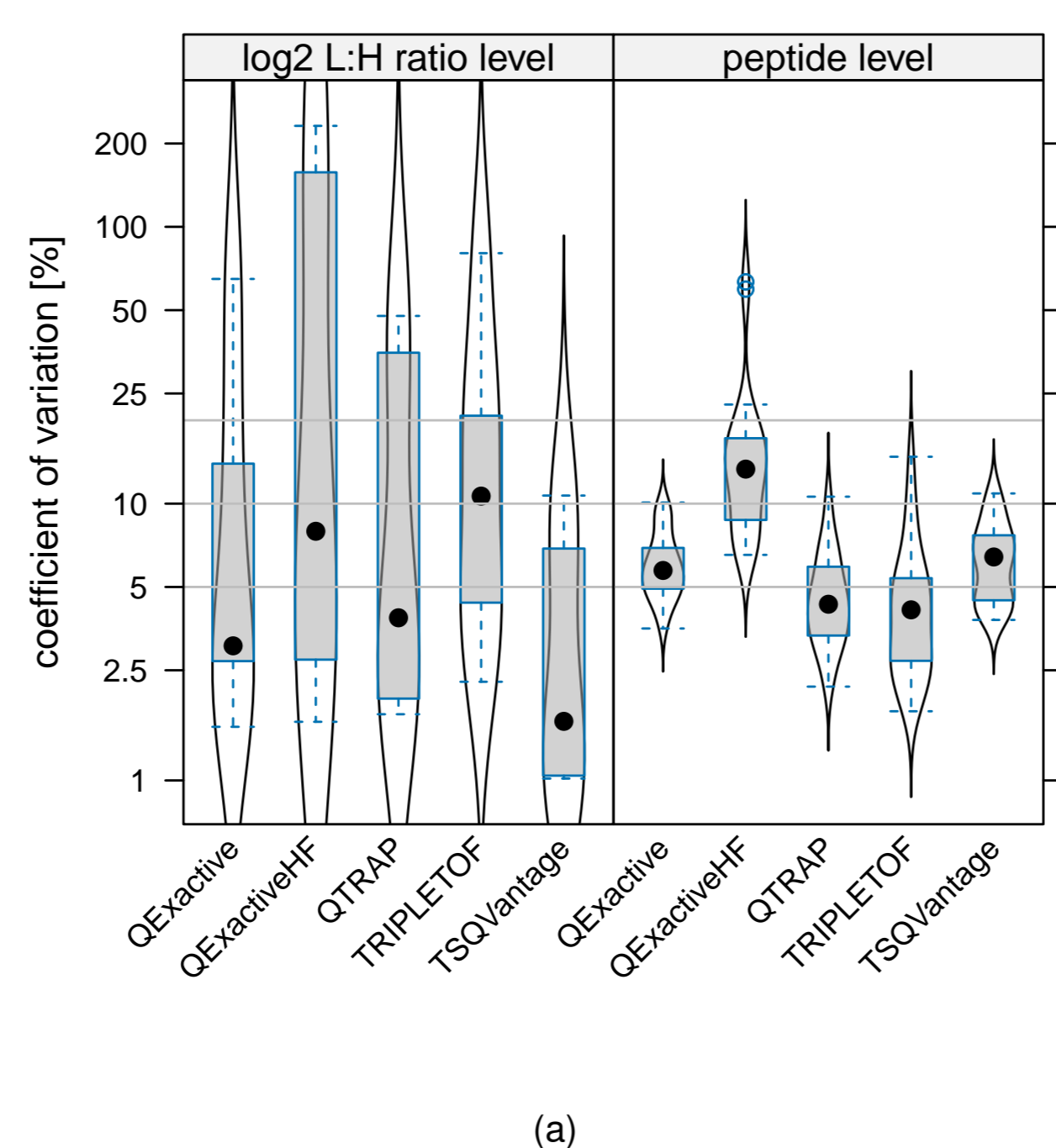


Figure 5: Precision – The violin and box plots display the distribution of the peptide level coefficient of variance (CV) computed on the light and heavy peptide signal of each platform.

## sigma mix peptide level signal

QExactive  
QExactiveHF  
QTRAP  
TRIPLETOF  
TSQVantage

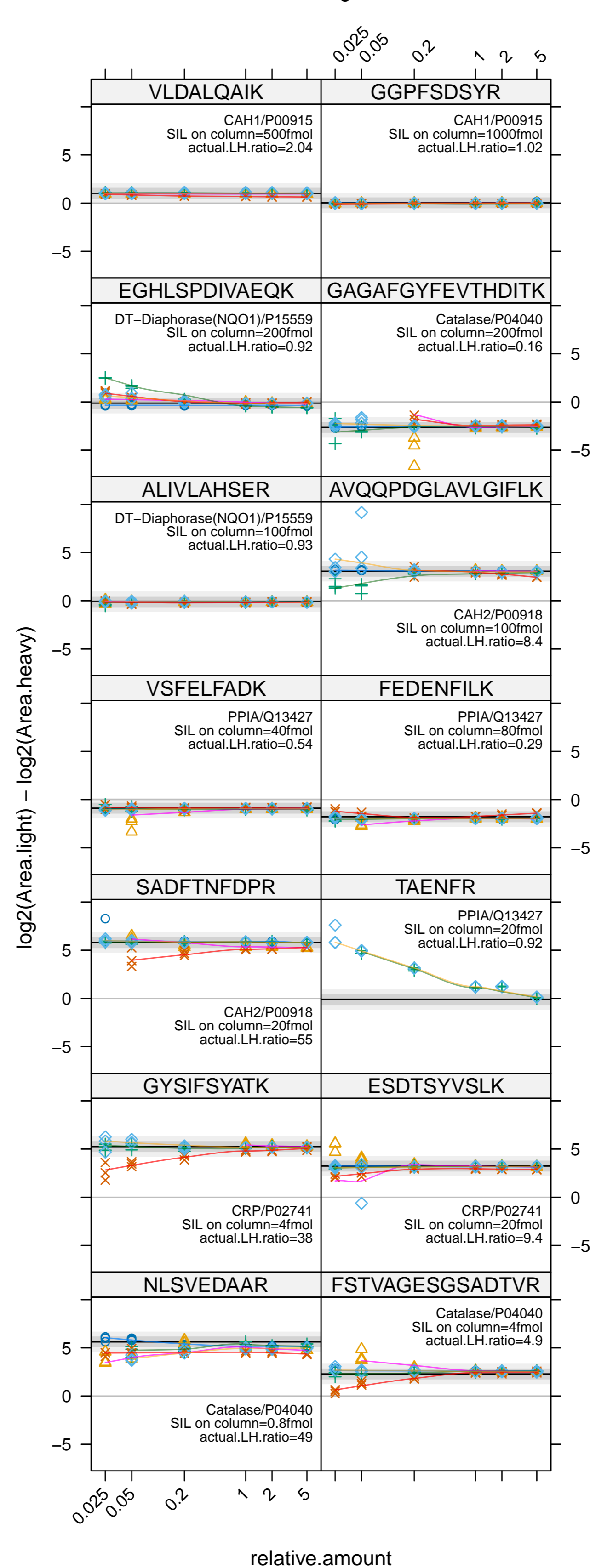

## log2 light heavy ratios of 6 dilutions on 5 MS platforms

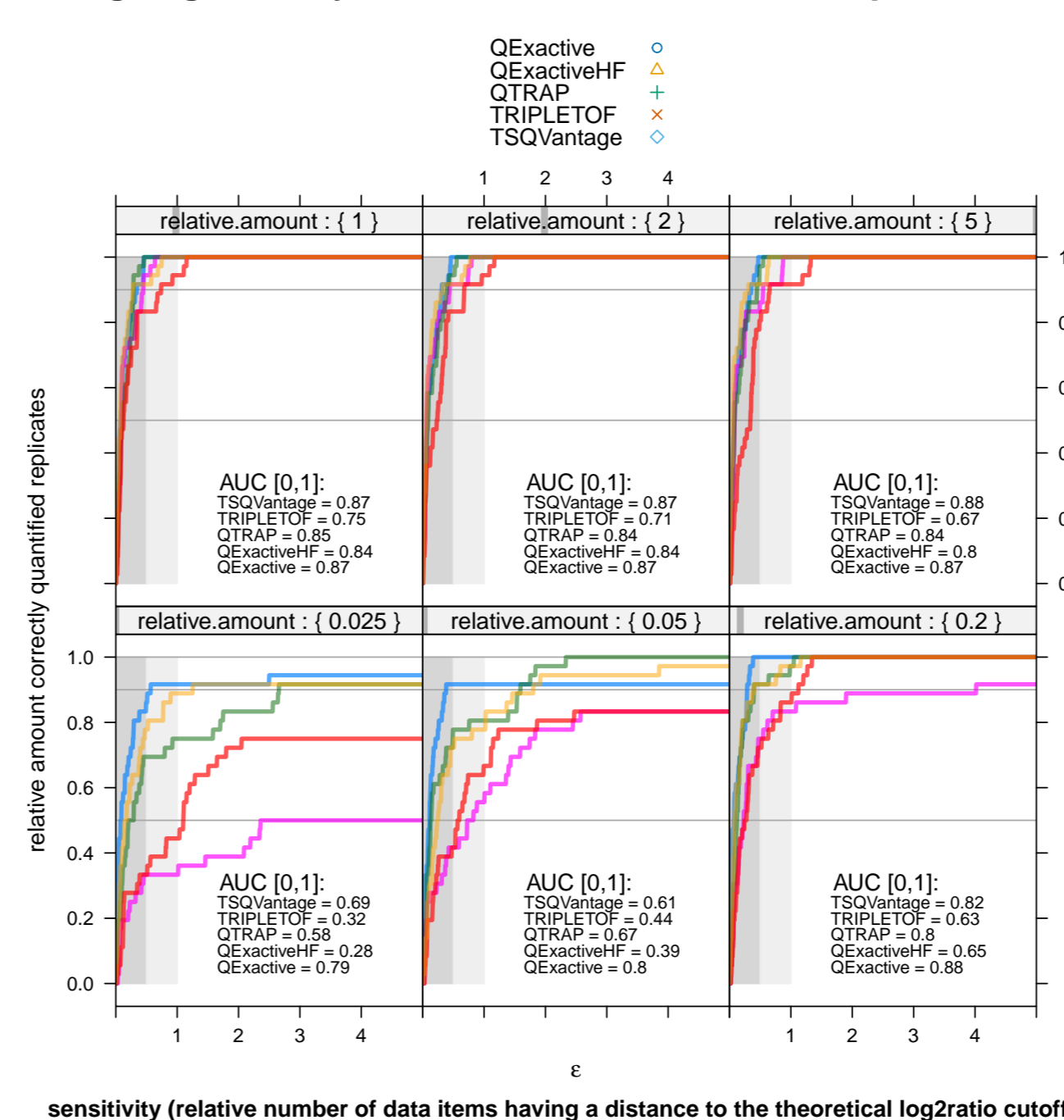


Figure 6: **Ratio stability upon analyte dilution and accuracy** – Upper chart: Each scatterplot panel displays the experimental derived  $\log_2$  L:H ratios versus the relative amount. Color grouping was done by instrument. The LOESS fit curves were added for visualizing the trend. The SIL value given in each panel legend is valid for the relative amount of 1. The horizontal black line indicates the theoretical  $\log_2$  L:H ratio. On both visualizations grey color boxes indicate the one and 2-linear-fold change. Lower chart: The graph displays in each panel a sensitivity curves for one relative amount for all the used workflows.



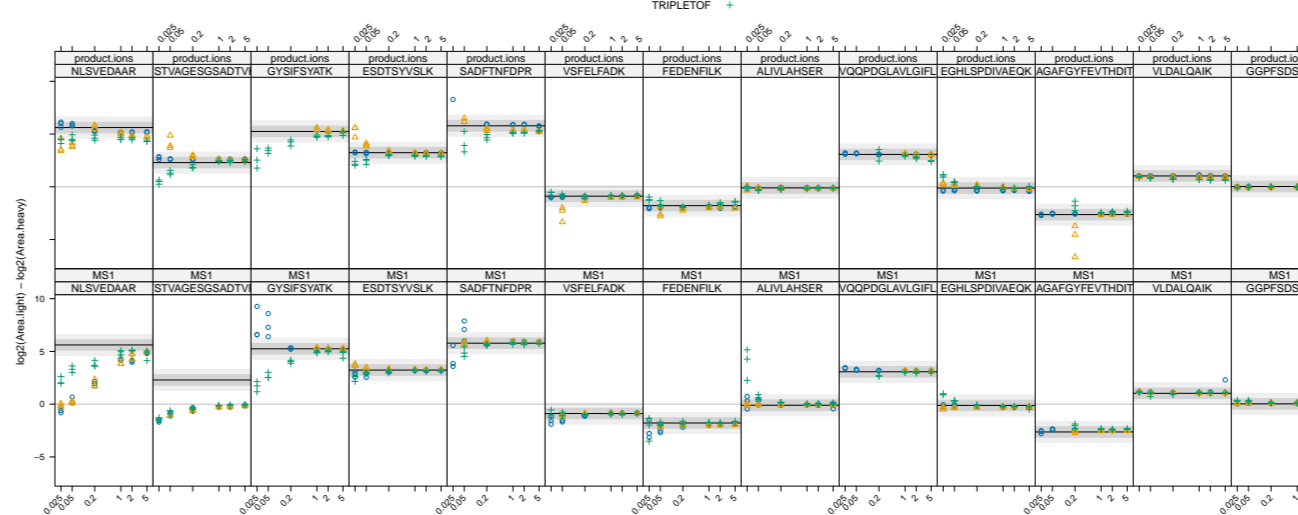


Figure 7: **Comparison of precursor and product ion** – Each panel displays the experimentally derived  $\log_2$  L:H ratios versus the relative amount. The black line indicates the reference  $\log_2$  L:H value. Dark and light grey shaded areas mark error margins of 1- and 2-fold on the linear scale, respectively.


## TSQ Vantage (SRM)

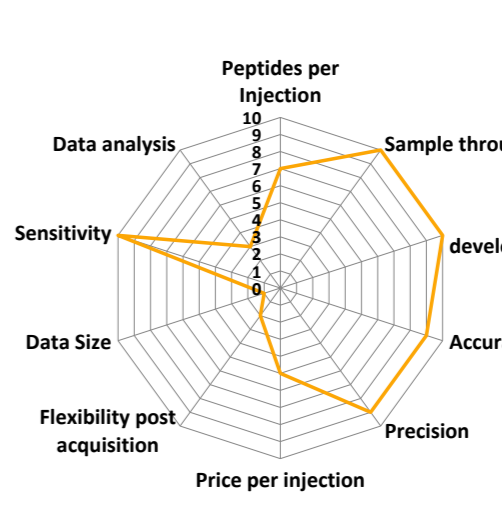


## QTRAP (SRM)

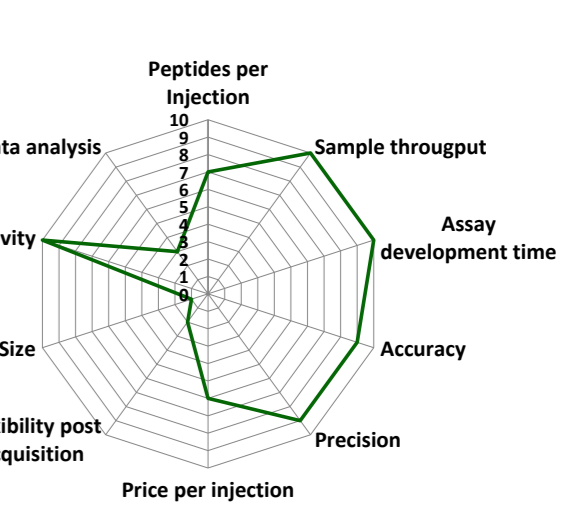

## QEXACTIVE (PRM)

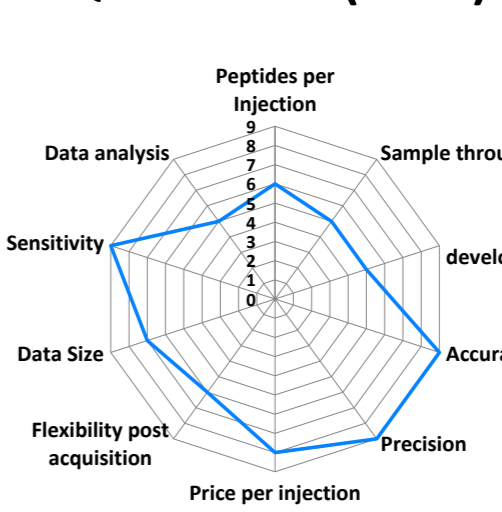

## QEXACTIVE HF (DIA)

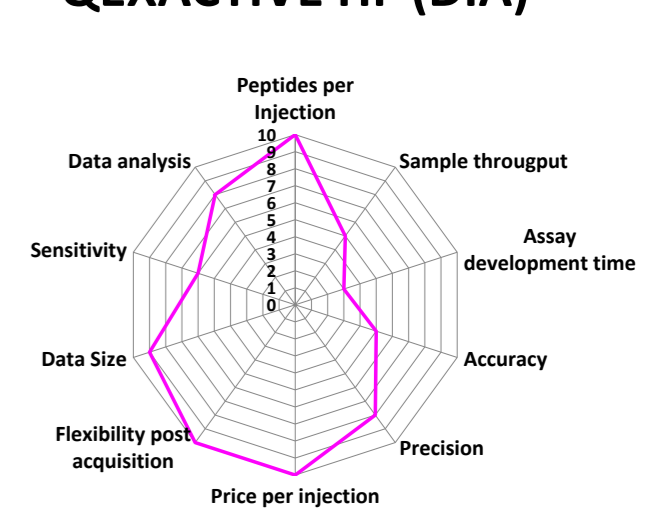

## TRIPLETOF (SWATH)

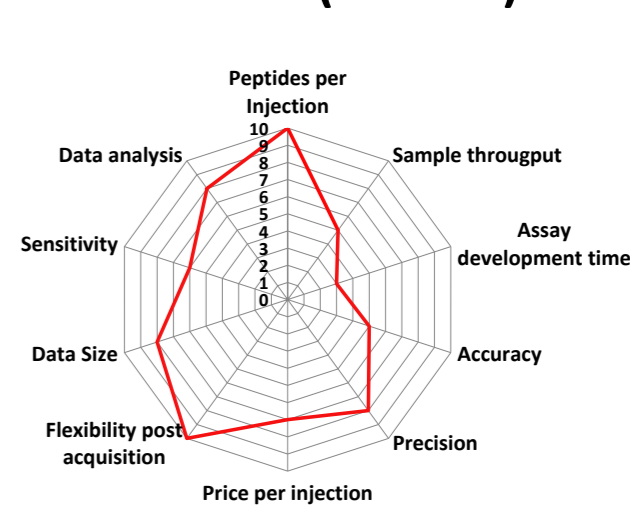


Figure 8: **Radarplot** of different workflow metrics scaled on an arbitrary axis with the units 0–10. Values do not necessarily represent global applicable values but rather reflect the situation at the FGZ. For the following categories the axis reflects 0 = worst, 10 = best: Peptides per injection, Sample throughput, Accuracy, Precision, Flexibility post acquisition, Sensitivity; 0 = best, 10 = worst: Assay development time, Price per injection, Data size, Data analysis.

## 4 Conclusions

The daily operational business of a core facility is normally not allowing for optimization of each method parameter, prior to recording a data set. This is reflected in this study by running the experiments with methods used on a routine basis in our core facility. Our study shows that targeted data acquisition (SRM, PRM) outperforms targeted data extraction strategies (DIA) with respect to quantitative accuracy and precision, especially when the analyte concentration are low. But targeted acquisition methods are only suited for monitoring a limited number of targets, while for screening experiments targeted data extraction workflows clearly outperform targeted data acquisition with respect to the number of features which can be followed. The question of whether greater throughput justifies lower sensitivity/specificity has yet to be answered in the context of the individual research project. In addition, large-scale DIA is very resource demanding (e.g. on the computational side) and especially data analysis needs more effort due to the high complexity of the raw data. Not all research environments are equally prepared for such challenges and should consider these aspects during experimental design setup.

**Acknowledgements** We thank Laura Kunz and Bernd Roschitzki for participating in our user study, and Bernd Roschitzki for critically reading our manuscript. We are grateful to Can Türker, Marco Schmidt and Ugur Gürel for providing the b-fabric information management system. We thank ETH Zurich and University of Zurich for their financial support.

## References

Kockmann T, Trachsel C, Panse C, Wähländer A, Selevsek N, Grossmann J, Wolski WE, Schlapbach R (2016). "Targeted proteomics coming of age - SRM, PRM and DIA performance evaluated from a core facility perspective." *PROTEOMICS*, pp. n/a–n/a. ISSN 1615-9861. doi:10.1002/pmic.201500502. URL <http://dx.doi.org/10.1002/pmic.201500502>.





